



CONTROL PANEL

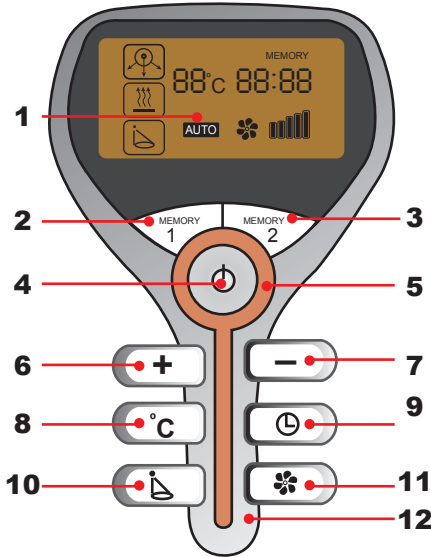

1. LCD Display Panel
2. MEMORY 1 Button: To set user's manual settings
(Initial setting: 104°F, 20 minutes)
3. MEMORY 2 Button: To set user's manual settings
(Initial setting: 111°F, 20 minutes)
4. Power Button: Power ON/OFF button
(Initial Setting: 108°F, 20 minutes, Fan: High, Laser: ON)
5. Power Lamp: Lights up when power is on.
6. Temperature/Time Up Button: Press button to increase output of far Infrared rays and raise temperature.
7. Temperature/Time Down Button: Press button to decrease output of Infrared rays and lower temperature.
8. Temperature Setting Button: Temperature can be set between 95°F~125°F.
9. Time Setting Button: Time can be set between 60 minutes.
10. Automatic Sensor Button: ON/OFF button for automatic movement
(When the Power button is on, this setting will not be indicated on the display panel.)
11. Fan: Fan takes in and takes out the air into the machine.
(Five settings(1~5) available for AUTO mode.)
Fan gets out stronger wind when you press the button .
12. Operation Indication Lamp: Indicates that the machine is running.

DISPLAY PANEL

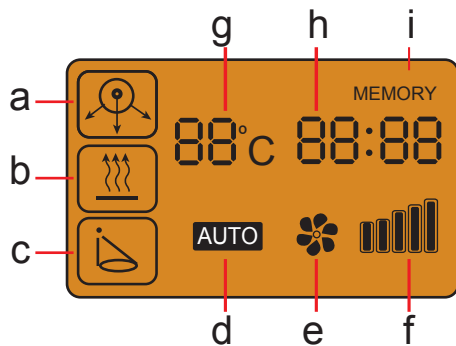

- a. Far Infrared Ray Indication
- b. Heat Therapy Operation Indication Icon
- c. Automatic Sensor Operation Indication Icon
- d. Automatic Fan Operation Indication Icon
- e. Fan Operation Indication Icon
- f. Fan Strength Indication Icon (5 bars: High, 3 bars: Medium, 1~5: Auto)
- g. Set Temperature Icon
- h. Set Time Icon
- i. Memory Setting Icon